

SIMPLE & EFFECTIVE FOR
LARGE STONE MANAGEMENT

EHL LITHOTRON EL27 AND ACCESSORIES

The Walz Lithotron EL27 Compact offers a secure, easy, and convenient way to administer bipolar ElectroHydraulic Lithotripsy (EHL).



KEY BENEFITS

CONTROL

Offering 3 different intensity options (up to 900 mJ), allowing users the ability to adjust energy levels based on stone densities.

SECURITY

Reduced tensile phases of concentrated energy enable accurate stone fragmentation while minimizing harm to soft tissues.


EFFICIENCY

The centered micro-shockwaves delivered by the probe facilitates targeting with fast pulses and sharply rising amplitudes.

EHL PROBES

- Highly flexible with enough tensile strength to maintain position.
- Precision machined probe tip for reproduceable shockwaves.
- Rounded probe tip edges for safe and smooth passage through endoscope.
- Centered micro spark plasmas preserve distances to the tissue and facilitates targeting.
- Multiple probe size options aiding in increased surface area contact with stone for efficient disintegration of concrements.

Electrohydraulic Lithotripsy (EHL)

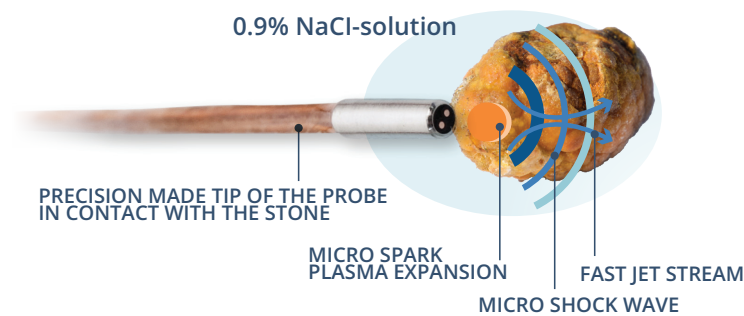

EHL PROBES

UPN/Order Number	REF/Catalog Number	Probe Diameter (Fr)	Probe Diameter (mm)	Working Length (cm)	Package Units
P243000FST	P243000FST	2.4	0.8	300	1/box
P303000FST	P303000FST	3.0	1.0	300	1/box
P453000FST	P453000FST	4.5	1.5	300	

*Note 3.0 Fr probe can be used with all 3.2 mm & 3.9 mm scope but recommended for 3.9 mm eyeMAX scope with 2.0 mm channel.

**Note 4.5 Fr probe to be used exclusively with 3.9 mm eyeMAX scope with 2.0 mm channel

EHL SYSTEM ACCESSORY

UPN/Order Number	REF/Catalog Number	Description	Package Units
FS2728B	FS2728B	Foot switch EL27 with connecting cable and connector	1/box
SKE180A	SKE180A	EHL Probe cable A2357 for EL27 Compact	1/box

CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.

