



TARGETED THERMAL HEMOSTASIS
**ENSURE
COAGULATION
FORCEPS**

EnSure™ reliable hemostasis of bleeding vessels without affecting surrounding tissue.

The EnSure Coagulation Forceps are designed to cauterize and coagulate GI bleeds. These single-use hemostatic forceps deliver targeted monopolar coagulation to isolate vessels and achieve hemostasis. Bleeding vessels can be approached with precision due to the reliable 360° bidirectional rotation.

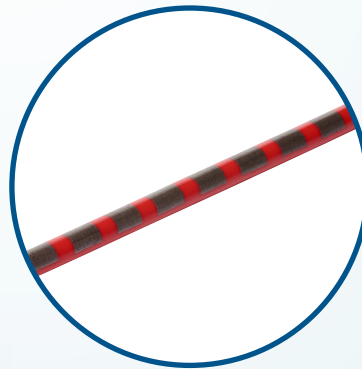

Markings on catheter to help determine insertion length within the endoscope.

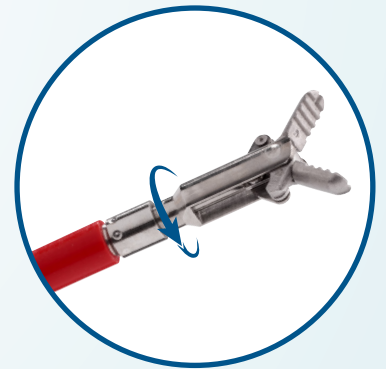


360° bidirectional rotation for targeting bleeds.

TWO JAW STYLES

Choose from a selection of two jaw styles to meet the clinical need:

- Tapered and Grooved
- Both types have jaws with teeth designed to prevent slippage during thermal coagulation.


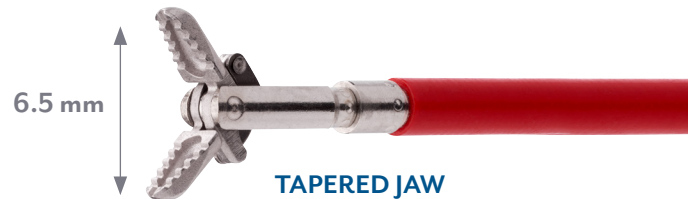


TAPERED JAW



GROOVED JAW

CROSS-SECTION OF A 2.8 MM WORKING CHANNEL

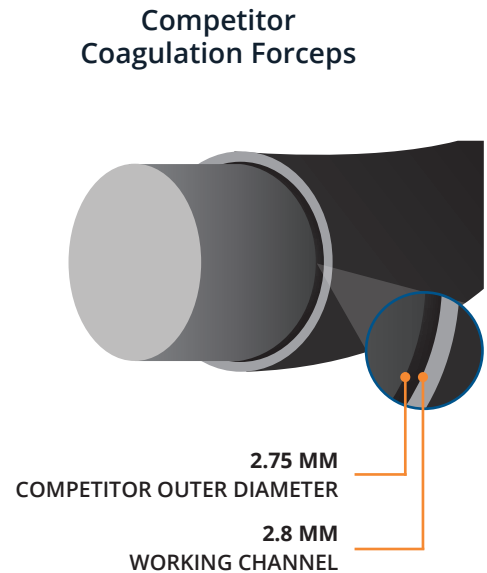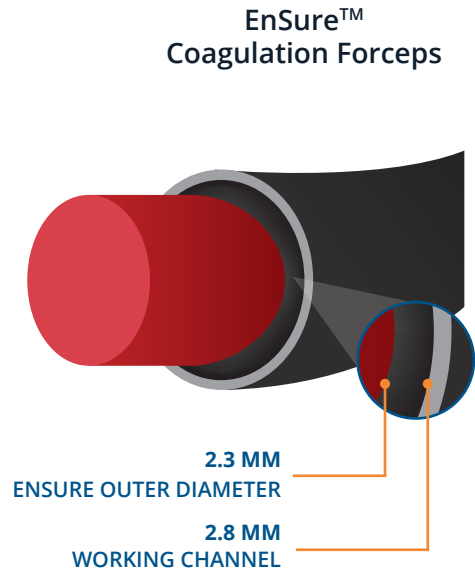

SLIM OUTER DIAMETER

The EnSure™ Coagulation Forceps feature a slim outer diameter of 2.3 mm that preserves approximately 0.5 mm of the working channel space for suction during hemostasis of bleeding vessels.

ENSURE COAGULATION FORCEPS STERILE

Order Number	Description	Jaw Style	Working Length (cm)	Outer Diameter (mm)	Jaw Opening (mm)	Package Units
CF48361	Ensure Coag Forceps Grooved Jaw	Grooved Jaw	230	2.3	6.3	10/Box
CF48331	Ensure Coag Forceps Grooved Jaw	Grooved Jaw	165	2.3	6.3	10/Box
CF48351	Ensure Coag Forceps Tapered Jaw	Tapered Jaw	230	2.3	6.5	10/box

Above graphical representation for illustrative purposes only.

CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.